

Accuracy calibration for BeneFusion VP5 infusion pump

1. Why should do accuracy calibration for infusion pump ?

The accuracy value defaulted as one brand pipe from manufacture .The accuracy parameter will be different due to different brand pipe in the end-users.

NOTES:


- 1).If you use your own brand of IV set, please calibrate the IV set on the machine before using.
- 2).After calibration, you can use same type of iv set on the machine.
- 3).It's better do calibration every unit with the IV Set.
- 4).If the pump installed drop counter, Please enter into the RATE MODE to adjust

2. How to do the accuracy calibration?

Step 1 Prepare the calibration tools.

- 1). Plug in AC Power (AC100-240V 50/60Hz)
- 2). Press POWER on (for 3-5 sec)
- 3). Screen automatically to enter the system
- 4).Prepare the measuring cylinder (5ml&50ml) and the Iv set that you need to adjust the accuracy.
- 5) One 50ml syringe needle

The view of calibration tools

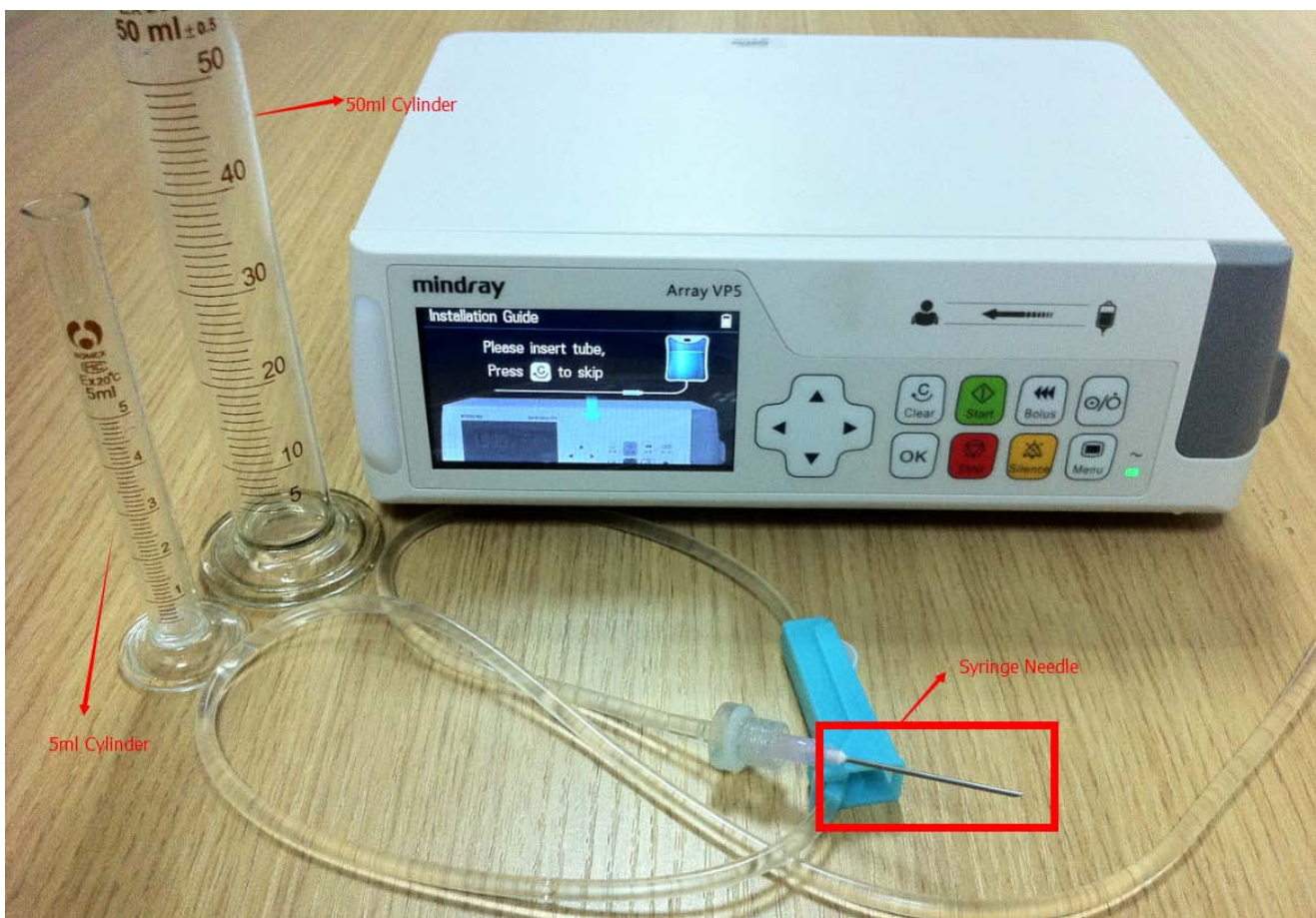

Step 2 Install the IV set



Pull up the door lock and install the IV set correctly (Keep the IV set straight)

Note: please remove the air bubble before accuracy calibration.





STEP 3 Enter the calibration interface and perform the calibration

STEP3-1 In the Main Menu, please choose the "System Maintenance" and input password "4321"

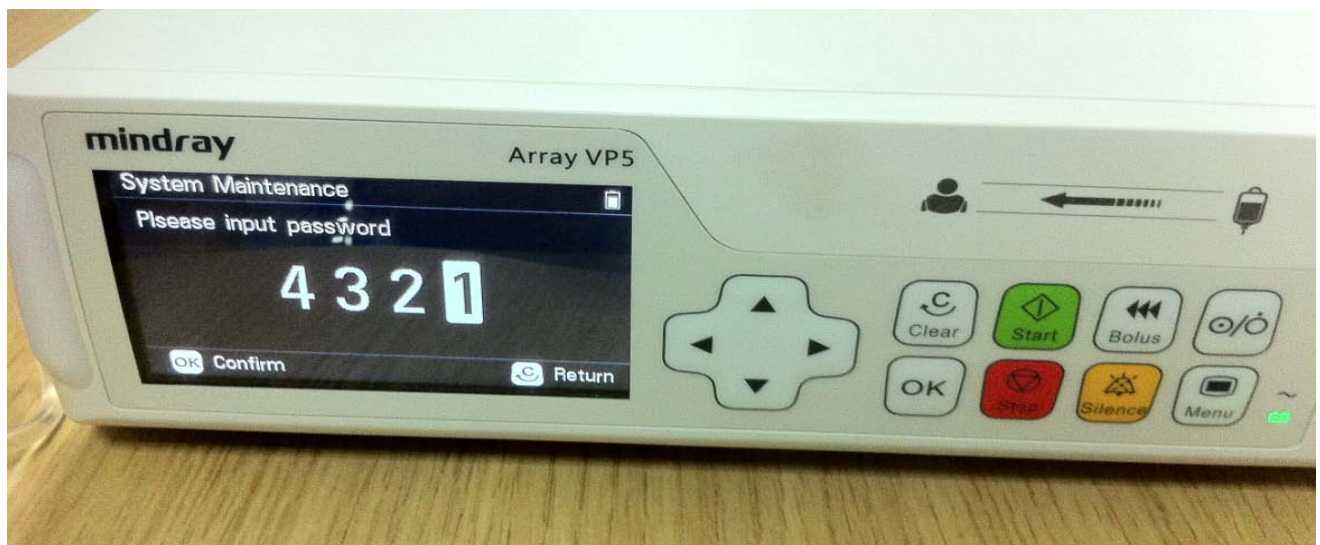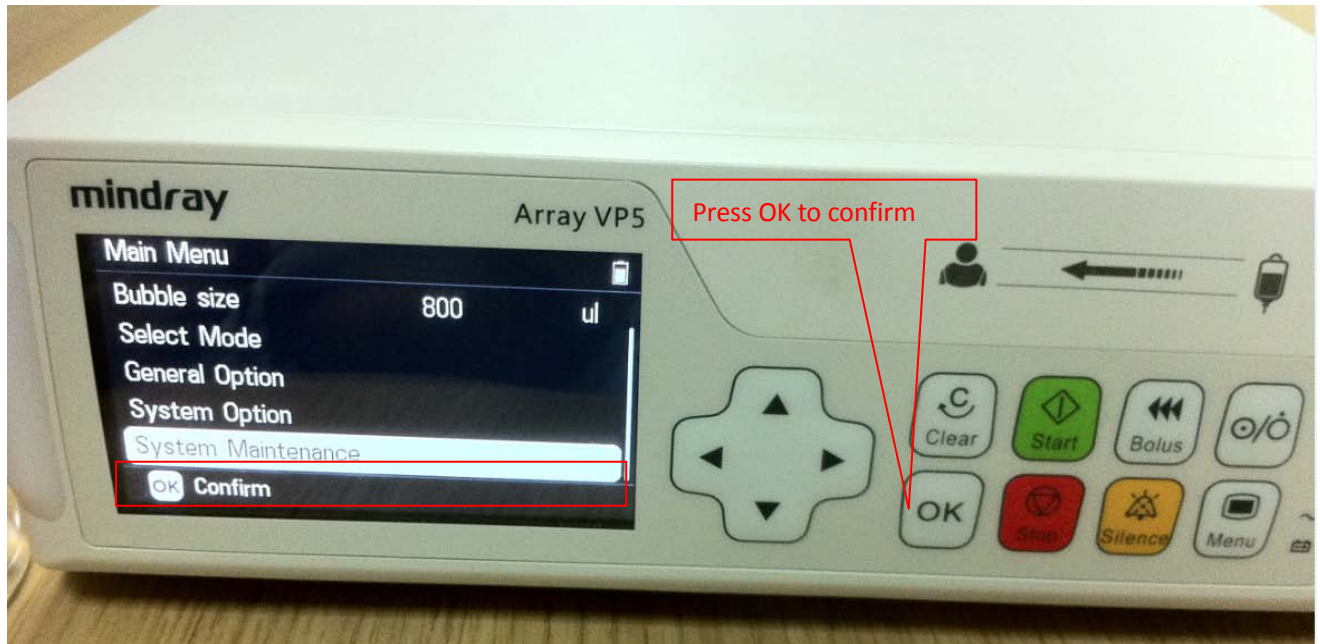



STEP 3-2 Choose the Option "Accuracy Calibration" and select the tube brand needs to be calibrated, for example "Brand 1"

Note: There are some brands of Iv set have been in-built in the pump, such as B.Broun, Shuangge. The accuracy already be calibrated, please choose to use this kinds of Iv set without calibration directly.

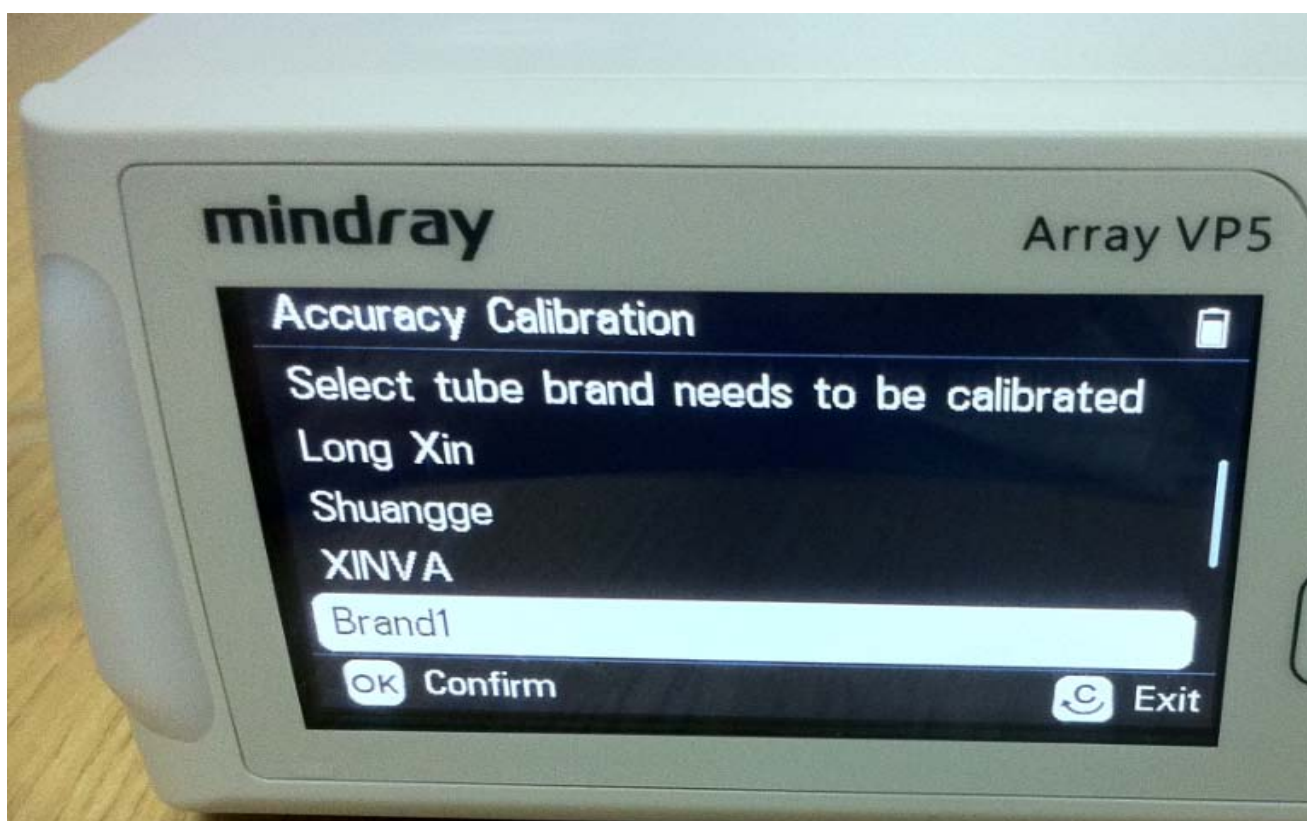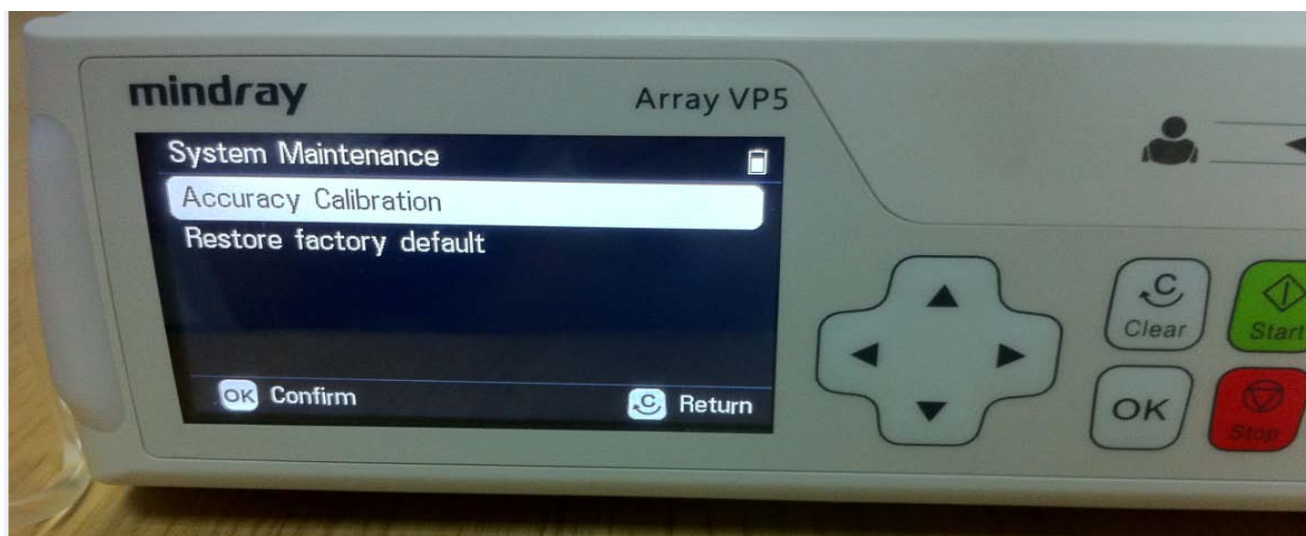



STEP 3-3 Finish the calibration procedure according to the instructions shows.

3-3-1 Use the [5 MI measuring cylinder](#) and recored the volume as below shown.

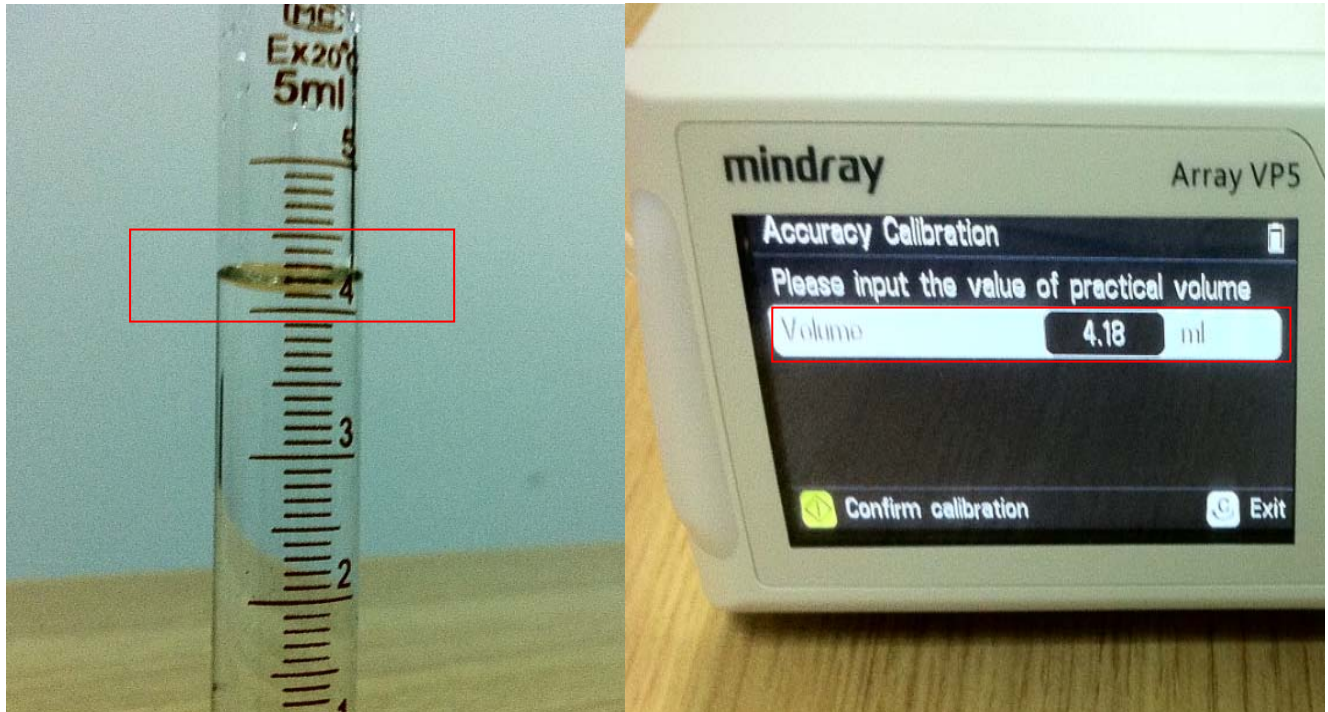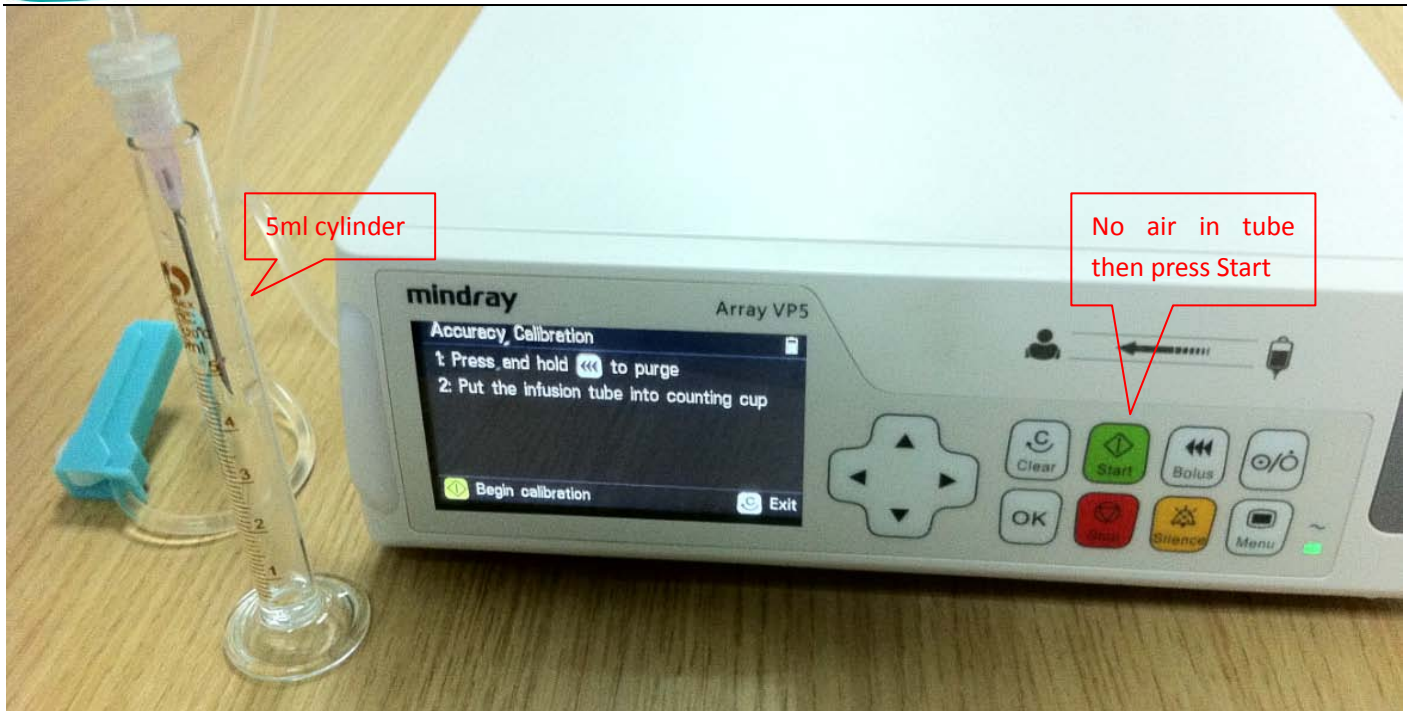

When it finish, input the liquid volume you read on the 5 ml measuring cylinder

3-3-2 Use the 50ml measuring cylinder and record the volume as below shown

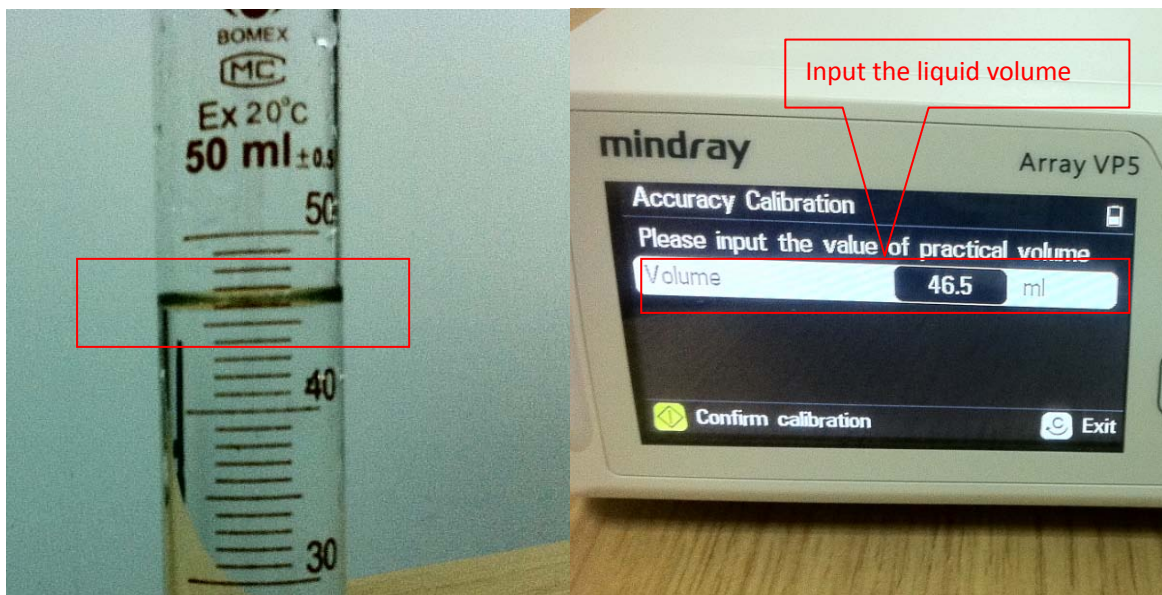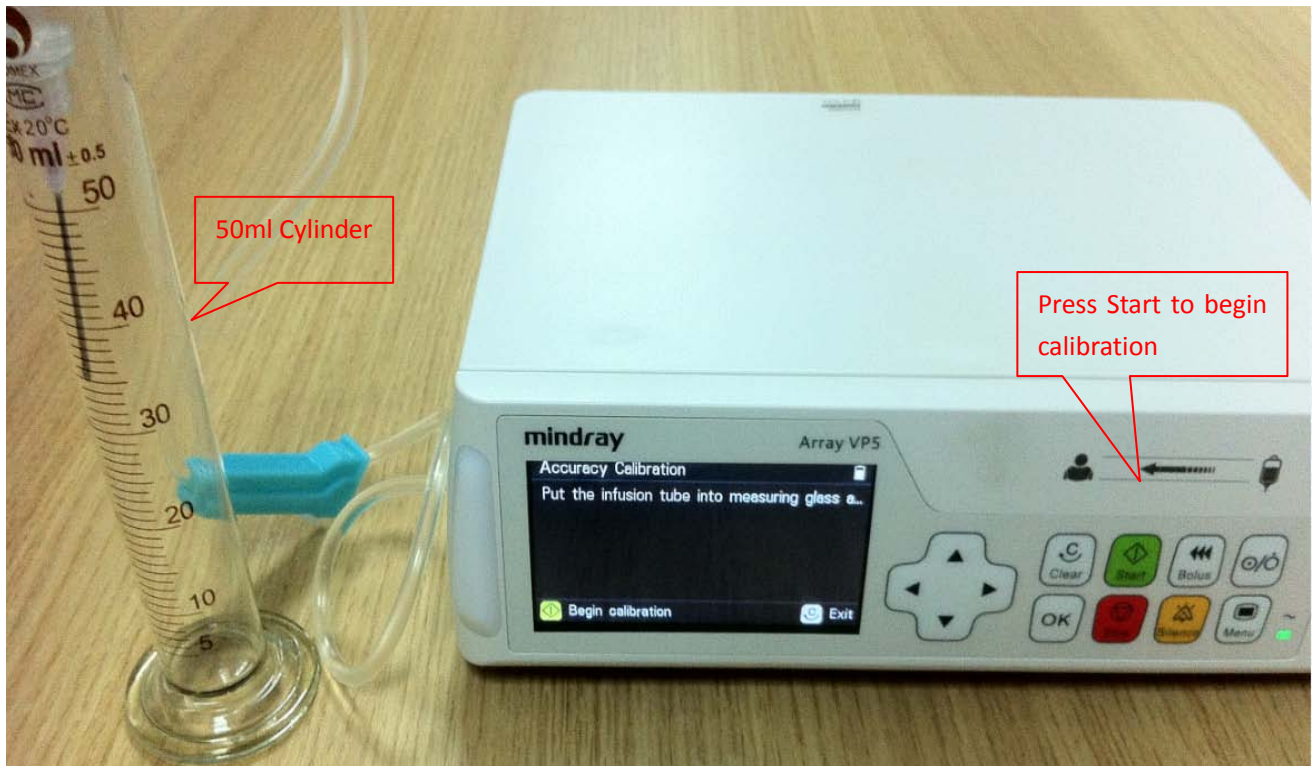

When it finish, input the liquid volume you read on the 50 ml measuring cylinder

STEP 4 After the steps above,we complete the accuracy calibration.



Press the Clear button to exit